
IN251200390V01_UK_FR_ES_DE_IT

700-199V70_700-199V90

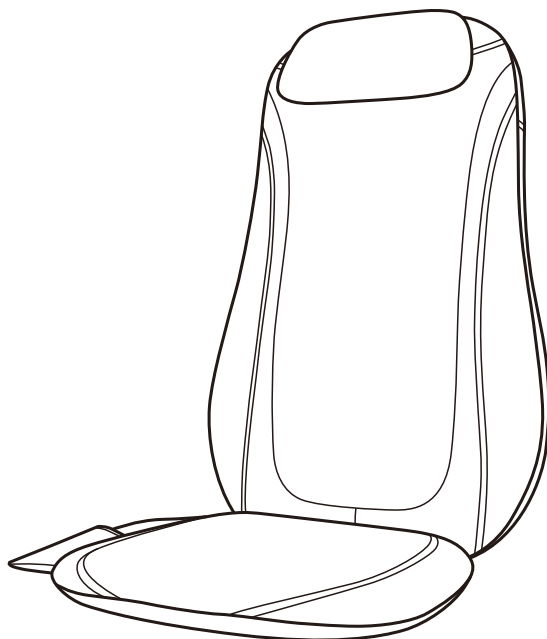Back Massager with Kneading and Heating

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EN_IMPORTANT, RETAIN FOR FUTURE REFERENCE: READ CAREFULLY.

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IMPORTANT SAFETY INSTRUCTIONS

RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use.

Explanation of symbols	
	This instruction manual belongs to this device. It contains important information about using this product. Read the instruction manual carefully. Non-observance of these instructions will result in serious injury or damage to the device.
	WARNING These warning notes must be observed to prevent any injury to the user.
	CAUTION These notes must be observed to prevent any damage to the device.
	NOTE These notes give you useful additional information on the installation or operation.

Important Safety Instruction

When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER- To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use and before cleaning.

WARNING

- To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not use and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Do not carry this product by its supply cord or use cord as a handle.
- 4) Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- 5) Do not drop or insert any object into any opening.
- 6) Do not use outdoors.
- 7) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- 8) To disconnect, switch off all controls and remove plug from outlet.
- 9) Close supervision is necessary when this appliance is used by, or near children, invalids, or disabled persons to ensure that they do not misuse the product.
- 10) Keep the cord away from heated surfaces. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation.
- 11)The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- 12)The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- 13) The appliance is only to be used with the power supply unit provided with the appliance.
- 14) This appliance is not intended for use by persons (including children)withreduced physical, Sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- 15)Special attention to components most susceptible to wear an advice to replace defective components immediately and/or keep the equipment out of use until repair.
- 16) Use this appliance only for its intended use as described in this manual.
- 17) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using this product.

Only use indoors!

Do not use the device in damp rooms (e.g. when bathing or showering)

- Do not wash!
- Do not dry clean!
- After using this product for 15 minutes, rest for 30 minutes, then use again.


SAVE THESE INSTRUCTIONS

INTRODUCTION

Thank you for purchasing the Back Massager with Kneading and Heating. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and Amazon order number.

INCLUDED

- 1 x Back Massager
- 1 x User Manual
- 1 x Car Charger
- 1 x Charger

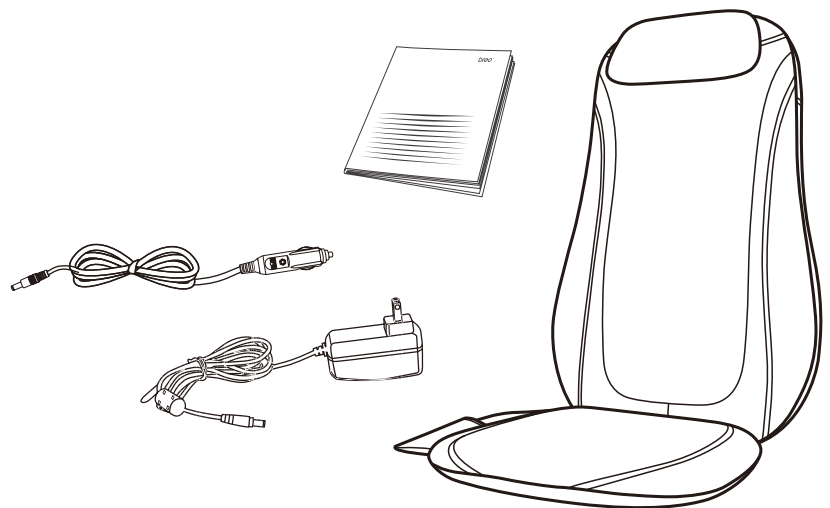
Product Specifications

Model	EMK-107G
Name	Back Massager
Power input	DC 12V == 2A
Power consumption	24 W
Default Time	15 mins
Enclosure material	PU, Nylon, Polyester, Latex, Resins
Output cable & Accessories	1.8m cable, Power Adaptor & Car Adapter

CE RoHS FC 

Made in China

NAME AND FUNCTION OF THE MACHINE PARTS

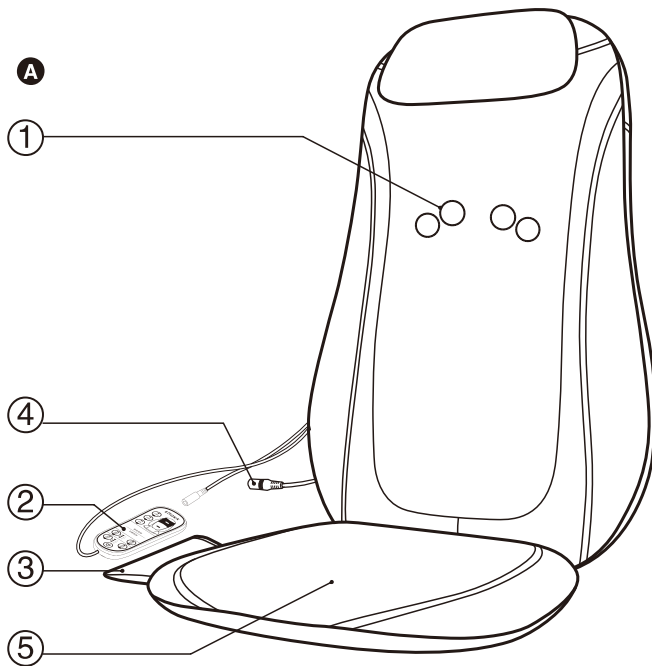
A. Main Unit

- ① Massage Head
- ② Remot Controller
- ③ Pocket for Remote Controller
- ④ DC connector
- ⑤ Seat Vibration

B. Power Adaptor

- ⑥ Plug
- ⑦ Cord
- ⑧ Appliance Plug

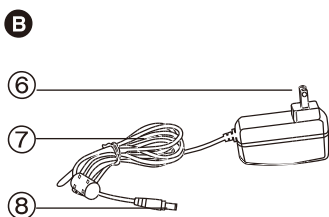
Vehicular cigar-lighter

HOW TO OPERATE

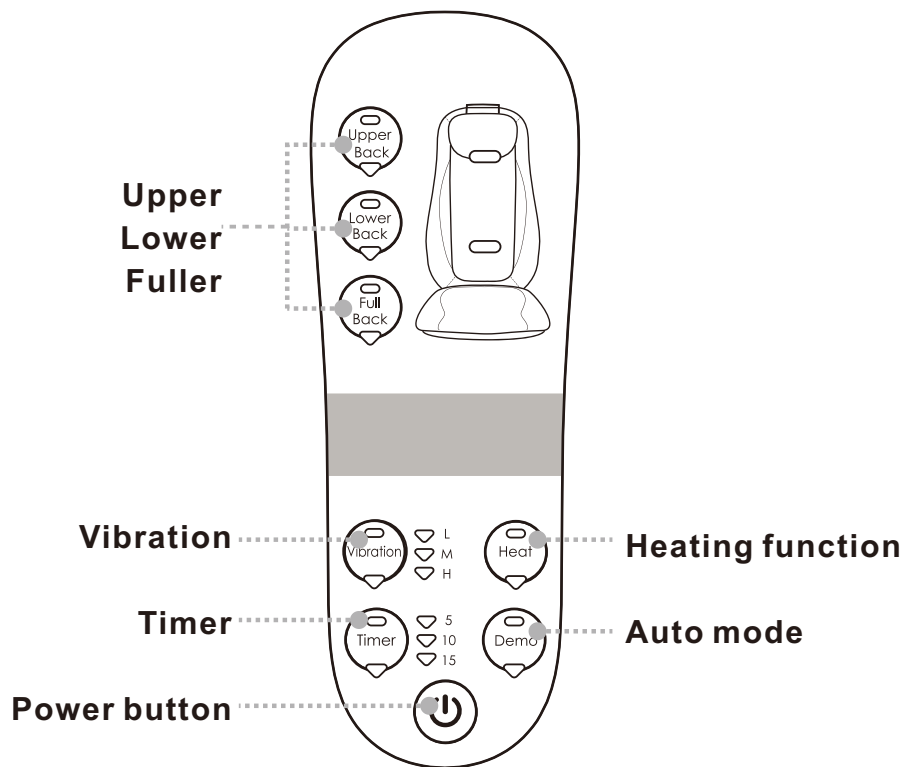Plug in the mains plug. Press the ON/OFF button, the appliance will begin to works. You can choose the button what you want, then the LED will light, the appliance will switch to the corresponding mode. If you want to stop the appliance, you can press the ON/OFF button. The massage head will reset to the initial position, waiting for the next job.

ATTENTION:

To ensure keep the appliance working properly, the DC connector must be fully inserted into the Appliance plug, until you hear the “click” sound.

	<p>ON/OFF Press the ON/OFF to start or stop the device.</p>
	<p>Demo Press DEMO button, it will demonstrates all the features and you can experience every feature. Pressing other button again will exit the demo mode, and enter normal working mode. If you press the DEMO button again in DEMO mode, the product will reset and the device will stop in 30S</p>
	<p>Timer Press timer button, set up the massage time for 5/10/15minutes.</p>

	<p>Upper Back While the device is working, press it to change the existing mode into Upper back mode.</p>
	<p>Lower Back While the device is working, press it to change the existing mode into Lower back mode.</p>
	<p>Full Back While the device is working, press it to change the existing mode into full back mode.</p>
	<p>Heat Press HEAT button, the product will be heating and the color of massage heads will be green & light orange. Press it again, the device will be stopped heating.</p>
	<p>Vibration Press vibration button, the seat will be vibrated</p>

TO CLEAN AND MAINTAIN THE PRODUCT

1. Cleaning

Before you clean, please make sure that the switch is in the off location, and the power plug has been pulled out from the socket. Please use dry cloth to clean the product.

2. Maintenance

Keep the product clean and far from the wet environment.

ATTENTION

Please do not use benzene or thinner and other aggressive detergents to clean the product or it may fade.


SUITABLE CROWD

Officer/Student
Using the computer for a long time leads to stiff muscles.

Elderly people
Muscle rigidity and backache.

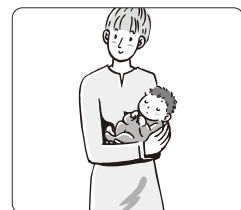
Housewife
Muscular soreness after busy housework.

Q&A

Problem	Cause	Remedy
Unable to start the device.	The cord is not connected properly.	Please push it harder until you heard the sound of "click".
The device doesn't work.	The usage time is too long.	Please stop using the product and unplug it from its power source for at least 1 hour.

Correct Disposal of this product

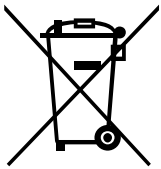This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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