

IN251200391V01_UK_FR_ES_DE_IT

700-200V00

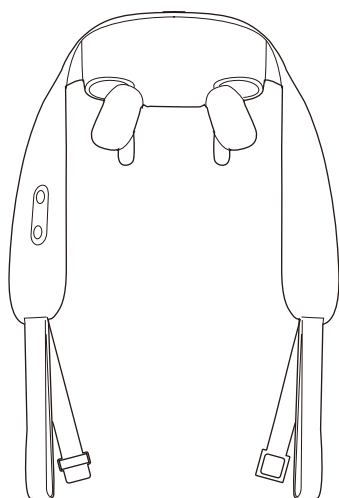3D kneading shoulder massager
Massager d'épaule à pétrissage 3D
Masajeador de hombros con amasado 3D
3D-Knet-Schultermassagegerät
Massaggiatore per spalle con impastamento 3D

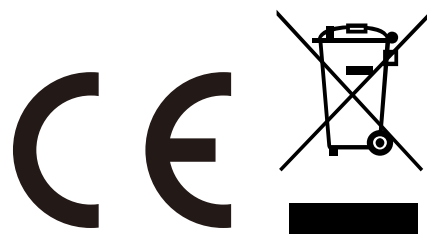EN_IMPORTANT, RETAIN FOR FUTURE REFERENCE: READ CAREFULLY.

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Please read this manual carefully before use and keep it in a safe place for future inquiry.

Important Safety Instruction

When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER- To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use and before cleaning.

WARNING

- To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not use and before putting on or taking of parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- 4) Do not drop or insert any object into any opening.
- 5) Do not use outdoors.
- 6) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- 7) Be sure to pull the plug out of the socket when you reach in. And don't pull it out when your hands are wet. That could cause electric shock and injury.
- 8) Do not carry this product by its supply cord or use cord as a handle.
- 9) This product contains batteries, it can only be replaced by technicians.
- 10) Remove the power plug from the socket immediately after the power failure. Otherwise, it may cause electric shock, short circuit and fire.
- 11) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for.
- 12) Special attention to components most susceptible to wear an advice to replace defective components immediately and/or keep the equipment out of use until repair.
- 13) Use this appliance only for its intended use as described in this manual.
- 14) Keep the cord away from heated surfaces. If the fault is found, stop using immediately and cut off the power supply. Avoid accidents or injuries. please do not repair it by yourself, please contact the seller or manufacturer in time, avoid accidents or injuries.
- 15) For charging use only protection class II power supply approved for use with household product, SELV output 5V, max. 2A.
- 16) The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.

- 17) This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- 18) Children should be supervised to ensure that they do not play with the appliance.
- 19) This appliance contains batteries that are non-replaceable.
- 20) The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

Please check whether the product and accessories are intact and damaged before using.

If you feel abnormal during use, please stop using immediately and consult your doctor. Otherwise may it be harmful to health.

1. People who are susceptible to electromagnetic malfunctions such as cardiac pacemakers and implant able medical electronic devices.
2. People who are supposed to lose their sense of temperature.
3. People with malignant tumors.
4. People with heart disease.
5. A person who is restless early in pregnancy or shortly after giving birth.
6. One who has an abnormal or curved spine.
7. A person with a disturbance of perception caused by a high degree of peripheral circulation in diabetes.
8. Someone with a scar on their skin.
9. People who need quiet and people who are in poor health.
10. People with osteoporosis, severe fractures, sprain, and muscle pain, etc., acute (head) pain.
11. People with 38 °C temperature (fever).
12. (acute inflammatory symptoms, tiredness, chills, changes in blood pressure, etc.) a time of intense, debilitating stress.
13. People with tendinitis.
14. Others treated at medical facilities.

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using this product.

Only use indoors!


Do not use the device in damp rooms (e.g. when bathing or showering)

- Do not wash!
- Do not dry clean!
- After using this product for 10 minutes, rest for 30 minutes, then use again.


SAVE THESE INSTRUCTIONS

Correct Disposal of this product


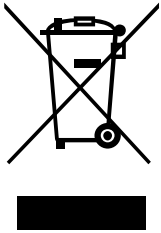
This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

NAME AND FUNCTION OF THE MACHINE PARTS

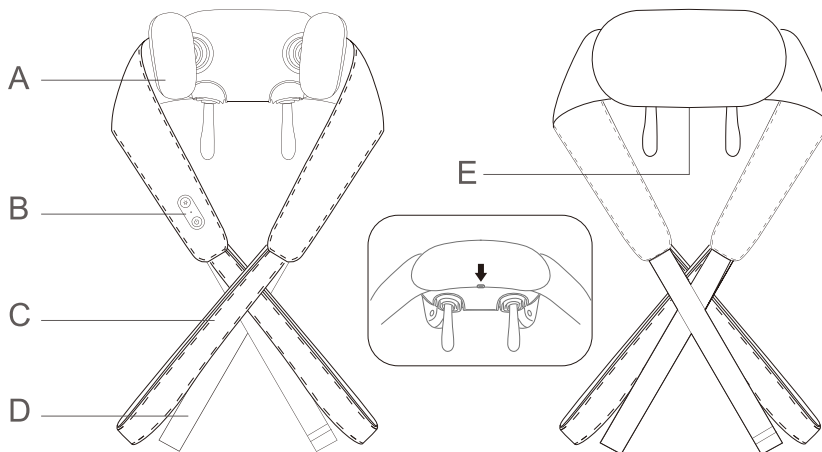- A. Massage node
- B. Control panel
- C. Hand Strap

- D. Elastic strap buckle
- E. Type-C charging port (on the back)

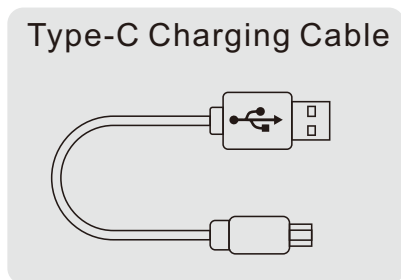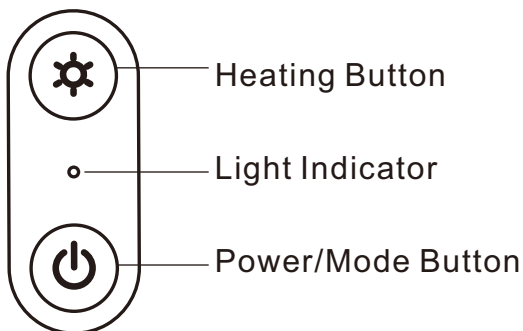Product Specifications


Model	EMK-158LQA
Name	3D kneading shoulder massager
Power input	DC 12V \approx 2A
Power consumption	18 W
Battery capacity	2000mAh


This product is suitable for home use only.

- When using this product must be equipped with a certified switch power cord to use.
- This product is only equipped with Type-C 5V charging data cable, if you need to charge the product, you need to match the charging head to charge, the charging head is recommended to use 5V2A specification charging head, the product built-in 7.4V/2000mAh lithium battery.
- If the product is not used for a long time, please charge the product within half a year to activate the battery.


MASSAGE POSITION

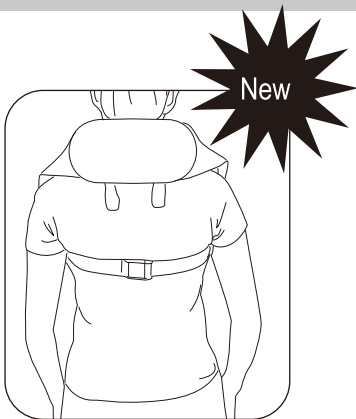**Neck
massage**

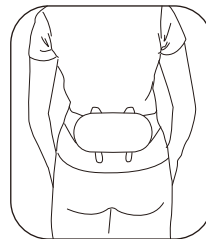**Lower back
relaxation**

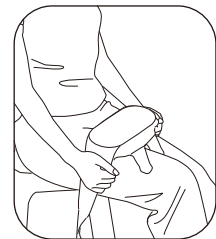**Relief leg
pressure**

Buckle strap design, completely free hands

Change the previous handle design, you can massage while working.

OPERATING INSTRUCTIONS

1. Power On: Long press the power button for 2 seconds to start the machine. The default mode is Kneading, indicated by a white light.
2. Mode Switching: Press the mode button to cycle through the two modes: Kneading Shiatsu Kneading.
3. Power Off: Long press the power button for 2 seconds to turn off the machine.
4. Automatic Shutdown: The machine will automatically shutdown after 10 minutes of continuous use
5. Heat Function
 - By default, the machine starts without heat.
 - Press the Heat Button to activate heat; the light turns red. Press again to turn off heat; the light returns to white.

HOW TO CHARGE

1. The machine cannot operate while charging.
2. When the machine is charging, the indicator light turns red.
3. When the battery is fully charged, the indicator light turns white.
4. If battery power drops below 20%, the indicator flashes red, signalling it's time to charge promptly.


Note: It is normal for the unit and charging plug to slightly heat up during charging.

Charging Time

The charging process takes approximately 3 hours.

After Charging

Once charging is complete, unplug the power plug from the outlet.

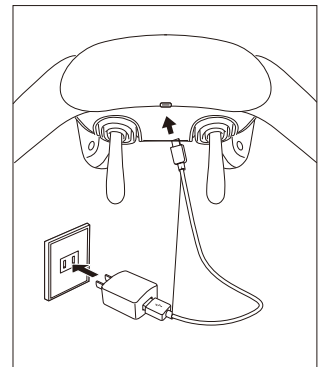TO CLEAN AND MAINTAIN THE PRODUCT

1. Cleaning

Before you clean, please make sure that the switch is in the off location, and the power plug has been pulled out from the socket.


Please use dry cloth to clean the product.

2. Maintenance

Keep the product clean and far from the wet environment.

Cleaning the main unit

- Wipe away dirt with a cloth soaked in soapy water, wringing out the cloth before wiping.
- Do not use alcohol, paint thinner, benzine, or nail polish remover. Using these may damage the appliance or cause cracks or discoloration in the components.

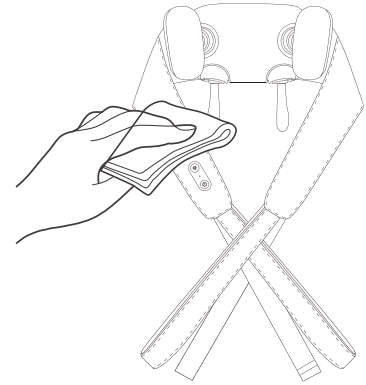OPERATIONS FOR USE

Problem	Cause	Remedy
Can not start.	The power button is not on.	Long press the power button to turn the massager on again.
	10 minutes automatic massage time is up.	Long press the power button to turn the massager on again.
	Excessive load placed on the product (e.g. someone standing on the product or heavy objects placed on the top of the product)	Remove excess load bearing items.
Sudden failure to start during use.	The battery of the product is exhausted.	Please charge in time.
There's a sound coming out of the product.	This is the sound of filling and discharging the internal air pump, which belongs to the normal working state.	no need to process.
The product is smoking.	The conductor of the product is short.	Discontinue use immediately and contact the designated service centre for this product.

If you do not use it for a long time, ensure that the function is normal before using it again.

periodical inspection	Are there any of the following faults	If the above fault occurs, in order to prevent accidents, please stop using immediately and cutoff the power supply. Send to designated maintenance point for repair.
	<ul style="list-style-type: none"> • The power cord is deeply bruised or deformed and smells burnt. • Cause abnormal sounds or vibrations during work. • After the switch is turned on, it does not work properly. • Leather case or appearance damaged. 	

If you still have any problem after using the methods above, please contact our customer support team.

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If you have any questions, please contact our customer care center, Our contact details are below:

US

 001-877-644-9366

 customerservice@aosom.com

Imported by Aosom LLC
27150 SW Kinsman Rd Wilsonville, OR 97070 USA
MADE IN CHINA

If you have any questions, please contact our customer care center, Our contact details are below:

CA

 416-792-6088

 customerservice@aosom.ca

Imported by Aosom Canada Inc.
7270 Woodbine Avenue, Unit 307, Markham,
Ontario Canada L3R 4B9
MADE IN CHINA

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If you have any questions, please contact our customer care center, Our contact details are below:

UK

 0044-800-240-4004

 enquiries@mhstar.co.uk

IMPORTER ADDRESS:
MH STAR UK Ltd
1 Northampton Cross Logistics Park
Northampton Cross Road
Northampton
NN4 9FH
MADE IN CHINA

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MH Handel GmbH