
IN251200393V01\_UK\_FR\_ES\_DE\_IT

**700-202V70\_700-202V90**

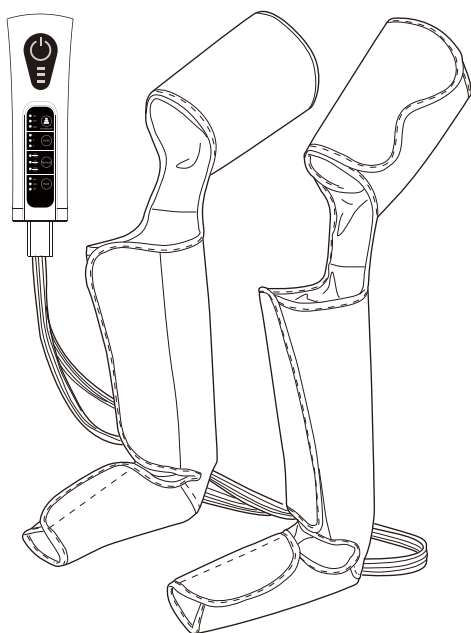**Air Compression Legs Massager**  
**Massageur de jambes par compression d'air**  
**Masajeador de piernas con compresión de aire**  
**Luftkompressions-Beinmassagegerät**  
**Massaggiatore per gambe a compressione d'aria**

**EN\_IMPORTANT, RETAIN FOR FUTURE REFERENCE: READ CAREFULLY.**

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## Important Safety Instruction

When using the product, all instructions and basic precautions should always be followed, including the following:

### **DANGER**

- To reduce the risk of electric shock: Always unplug this product from the electrical outlet immediately after use and before cleaning.

### **WARNING**

- To reduce the risk of burns, fire, electric shock, or injury to persons:
  - 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not use and before putting on or taking off parts.
  - 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
  - 3) Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
  - 4) Do not drop or insert any object into any opening.
  - 5) Do not use outdoors.
  - 6) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
  - 7) Be sure to pull the plug out of the socket when you reach in. And don't pull it out when your hands are wet. That could cause electric shock and injury.
  - 8) Do not carry this product by its supply cord or use cord as a handle.
  - 9) Keep the cord away from heated surfaces. If the fault is found, stop using immediately and cut off the power supply. Avoid accidents or injuries. please do not repair it by yourself, please contact the seller or manufacturer in time, avoid accidents or injuries.
  - 10) Remove the power plug from the socket immediately after the power failure. Otherwise, it may cause electric shock, short circuit and fire.
  - 11) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination.
  - 12) Special attention to components most susceptible to wear an advice to replace defective components immediately and/or keep the equipment out of use until repair.
  - 13) Use this appliance only for its intended use as described in this manual.
  - 14) The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
  - 15) The appliance is only to be used with the power supply unit provided with the appliance.
  - 16) As the surface of this product may get hot, people who are sensitive to heat should pay attention to using it to avoid minor burns.

17) This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

**If you feel abnormal during use, please stop using immediately and consult your doctor. Otherwise may it be harmful to health.**

1. People who are susceptible to electromagnetic malfunctions such as cardiac pacemakers and implant able medical electronic devices.
2. People who are supposed to lose their sense of temperature.
3. People with malignant tumors.
4. People with heart disease.
5. A person who is restless early in pregnancy or shortly after giving birth.
6. One who has an abnormal or curved spine.
7. A person with a disturbance of perception caused by a high degree of peripheral circulation in diabetes.
8. Someone with a scar on their skin.
9. People who need quiet and people who are in poor health.
10. People with osteoporosis, severe fractures, sprain, and muscle pain, etc., acute (head) pain.
11. People with 38 °C temperature (fever).
12. (acute inflammatory symptoms, tiredness, chills, changes in blood pressure, etc.) a time of intense, debilitating stress.
13. People with tendinitis.
14. Others treated at medical facilities.

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using this product.

Only use indoors!

Do not use the device in damp rooms (e.g. when bathing or showering)

Do not wash!

Do not dry clean!

**SAVE THESE INSTRUCTIONS**

## Introduction

Thank you for purchasing the Air Compression Massager for Feet, Leg and Thigh. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number.

## Package Includes

2 x Leg Wrap ("R" for right leg/"L" for left leg)

1 x Adapter

1 x User Manual

1 x Handheld Controller

## Product Specifications


Model	EMK-703D
Name	Air Compression Legs Massager
Power input	DC 12V $\approx$ 2A
Default Time	10 mins

This product is intended for house-hold use only.

Short time operation 10 minutes, resting time 30 minutes.


**Note:** Certified manufacturer's adapter must be equipped in the use of this products.

## Product Diagram

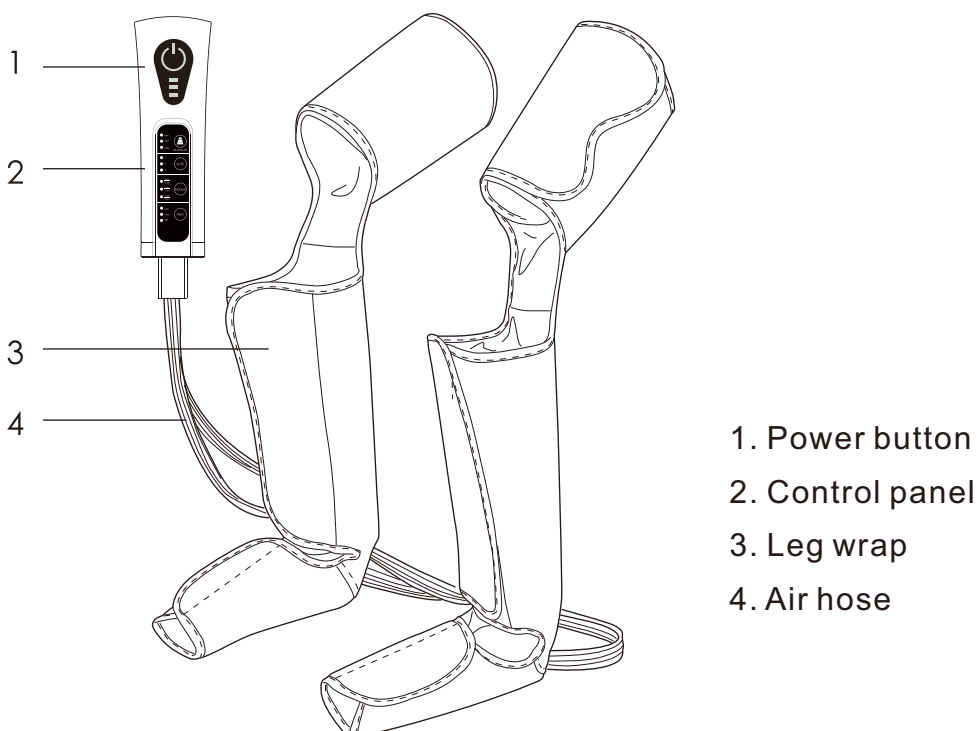
1. Power button
2. Control panel
3. Leg wrap
4. Air hose


## Accessories

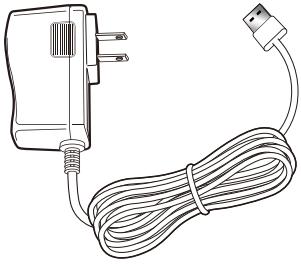Power Adapter

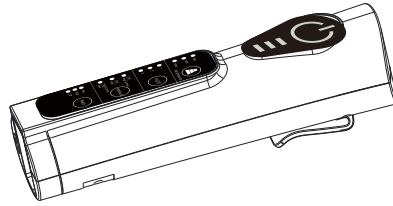Remote Controller

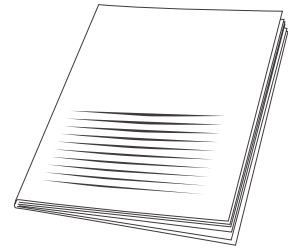
User Manual


## Product Operation Guide

### Before Using


1. Inset the air nozzle into the 7 corresponding holes of remote controller.

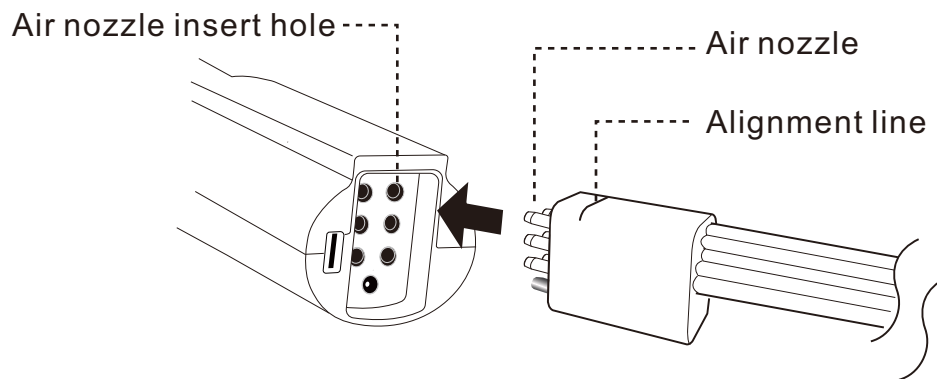
2. Insert the USB interfaces of adapter into remote controller's power hole, the other end of adapter connect with power supply.

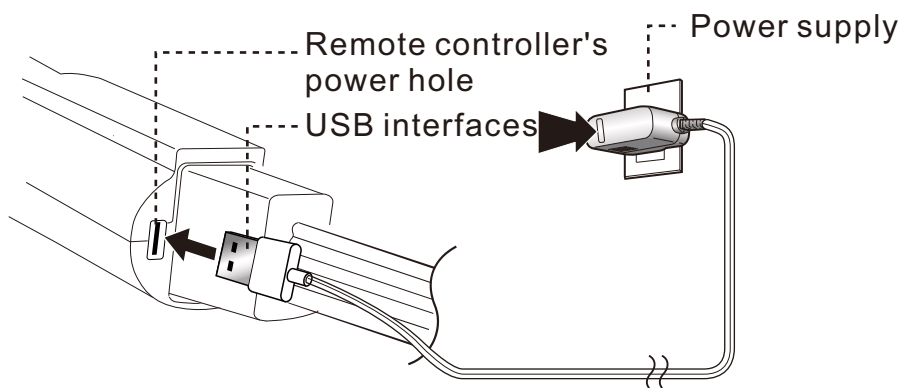1. Flat your feet and put it into the correct cover. "R" is the right foot, and "L" the left foot.
2. Paste the foot part tightly.
3. Paste the leg & thigh part tightly.
4. Move your feet and make sure your feet are tied in the right place.
5. Insert the air pressure connectors into the corresponding holes on controller.

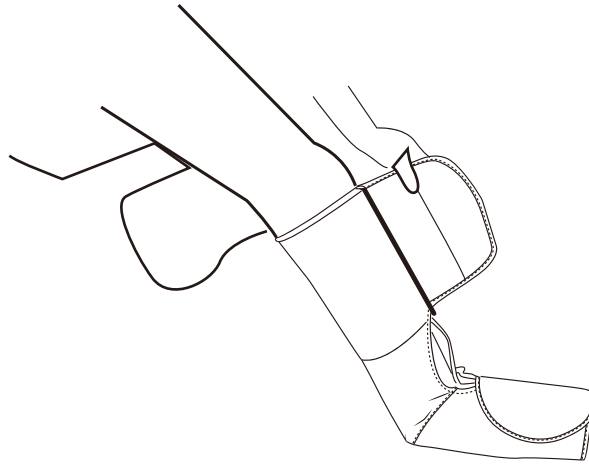Feet, Leg & Thigh

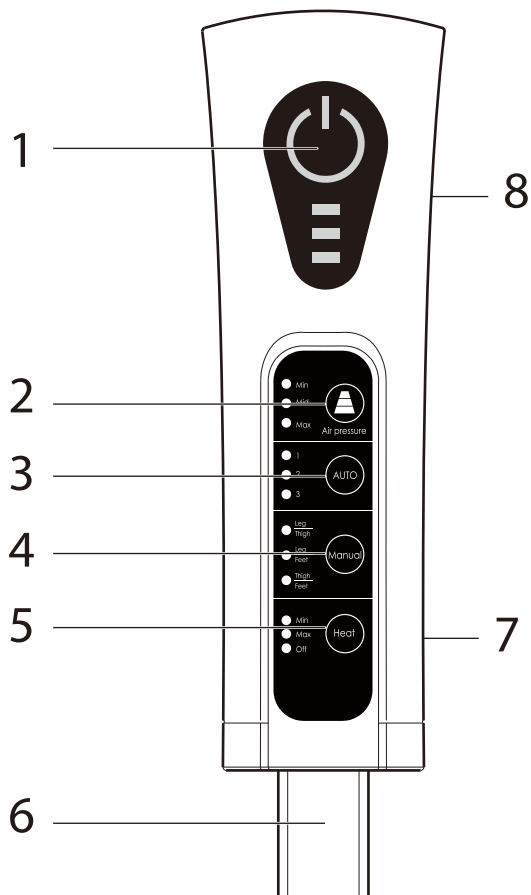1. Power button
2. Air pressure intensity button
3. Auto button
4. Manual button
5. Heating button
6. Air channel
7. Bottom shell
8. Hanging hook

**Power button:**

Press the Power button to start the product, default mode is Auto 1 and Mid intensity. Press the Power button again to turn off the product.

**Air pressure intensity button:**

When the product is working, press the button to change the intensity as: Max level→Min level→Mid level→Max level, in turn. The indicator light will be on when its changing.

**Auto button:**

When the product is working, press the Auto button to change the modes as: Auto2 →Auto3→Auto1→Auto2, in turn. The indicator light will be on when its changing.

**Manual button:**

When the product is working, press the Manual button to start the Manual mode as: (Leg & Thigh)→(Leg & Feet)→(Thigh & Feet)→ (Leg & Feet), in turn. The indicator light will be on when its changing.

**Heating button:**

When the product is working, heating function is not working and the indicator Light(Off) lights by default.


Press the button to start the heating function as: Min level→Max level→Off heat→Min level, in turn. The indicator light will be on when its changing.

**Note:** The default working time of the product is 10 minutes; If the product has worked continuously for 10 minutes, all massage movements would stop and then it would shut off automatically.

**Q&A**

Problem	Cause	Remedy
The product could not be turned on normally.	The power button has not been pressed.	Please press the power button.
The product gets hot.	The product has been used continuously for long.	Please wait for 30 minutes before using again.
There exists air leakage.	The air hoses are not installed properly.	Please install the air hoses properly.

## Correct Disposal of this product

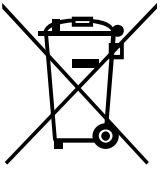This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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