

IN251200394V01_UK_FR_ES_DE_IT

700-203V70_700-203V90

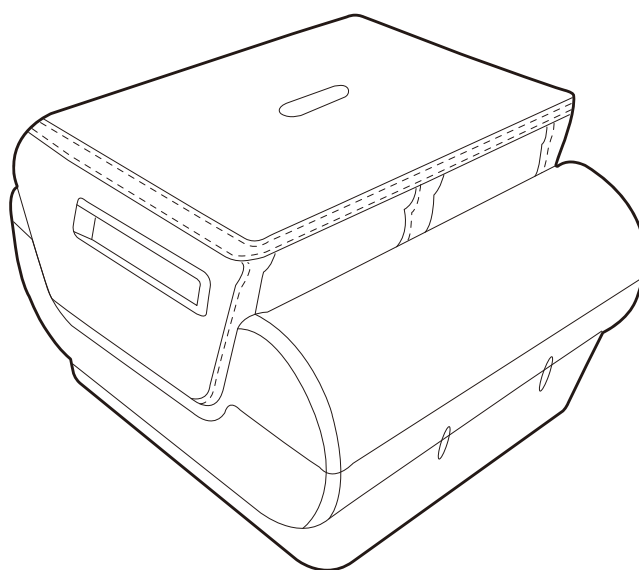Foot Massager
Appareil de massage des pieds
Masajeador de pies
Fußmassagegerät
Massaggiatore piedi

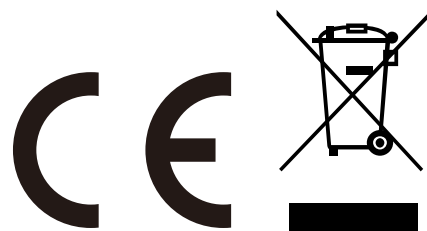EN_IMPORTANT, RETAIN FOR FUTURE REFERENCE: READ CAREFULLY.

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Important Safety Instruction

When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER

- To reduce the risk of electric shock: Always unplug this product from the electrical outlet immediately after use and before cleaning.

WARNING

- To reduce the risk of burns, fire, electric shock, or injury to persons:
 - 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not use and before putting on or taking off parts.
 - 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
 - 3) Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
 - 4) Do not drop or insert any object into any opening.
 - 5) Do not use outdoors.
 - 6) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
 - 7) Be sure to pull the plug out of the socket when you reach in. And don't pull it out when your hands are wet. That could cause electric shock and injury.
 - 8) Do not carry this product by its supply cord or use cord as a handle.
 - 9) Keep the cord away from heated surfaces. If the fault is found, stop using immediately and cut off the power supply. Avoid accidents or injuries. please do not repair it by yourself, please contact the seller or manufacturer in time, avoid accidents or injuries.
 - 10) Remove the power plug from the socket immediately after the power failure. Otherwise, it may cause electric shock, short circuit and fire.
 - 11) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination.
 - 12) Special attention to components most susceptible to wear an advice to replace defective components immediately and/or keep the equipment out of use until repair.
 - 13) Use this appliance only for its intended use as described in this manual.
 - 14) The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
 - 15) The appliance is only to be used with the power supply unit provided with the appliance.
 - 16) As the surface of this product may get hot, people who are sensitive to heat should pay attention to using it to avoid minor burns.

17) This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

If you feel abnormal during use, please stop using immediately and consult your doctor. Otherwise may it be harmful to health.

1. People who are susceptible to electromagnetic malfunctions such as cardiac pacemakers and implant able medical electronic devices.
2. People who are supposed to lose their sense of temperature.
3. People with malignant tumors.
4. People with heart disease.
5. A person who is restless early in pregnancy or shortly after giving birth.
6. One who has an abnormal or curved spine.
7. A person with a disturbance of perception caused by a high degree of peripheral circulation in diabetes.
8. Someone with a scar on their skin.
9. People who need quiet and people who are in poor health.
10. People with osteoporosis, severe fractures, sprain, and muscle pain, etc., acute (head) pain.
11. People with 38 °C temperature (fever).
12. (acute inflammatory symptoms, tiredness, chills, changes in blood pressure, etc.) a time of intense, debilitating stress.
13. People with tendinitis.
14. Others treated at medical facilities.

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using this product.

Only use indoors!

Do not use the device in damp rooms (e.g. when bathing or showering)

Do not wash!


Do not dry clean!

SAVE THESE INSTRUCTIONS


Product structure

Product body


- 1. Control panel
- 2. Removable and washable cloth cover
- 3. Main body
- 4. Non-slip foot pad
- 5. Foldable

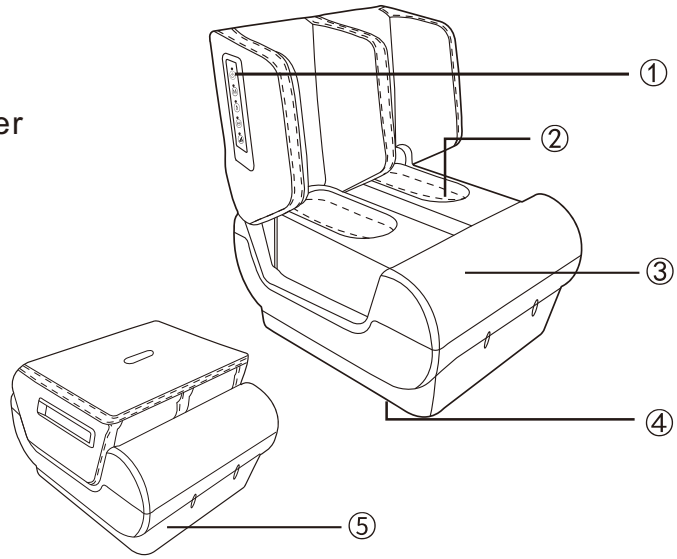Product bottom

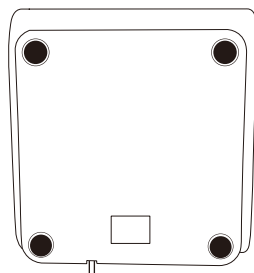Power usage diagram

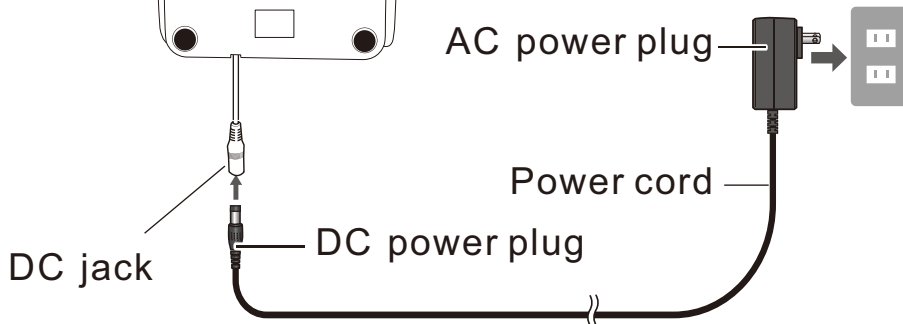ATTENTION

When using this product, be sure to use it with our company's certified switching power supply.

Instructions

Power button

Mode button

Heating button

Air pressure switch

Air pressure intensity

Connect the DC jack on the main unit to the DC power plug. When the AC plug is powered on, the product will emit a beep-ing sound.

Power button (function start/stop)

- In the off state, press this button to turn on the power, and the corresponding green LED will light up.
- In the on state, press this button to turn off the power, and the corresponding green LED will go out.
- After turning on, the default timer is 15 minutes, and it auto-matically enters the off state.
- There is no output after turning on, and the mode button, air pressure button or heating button must be operated to respond.

Air pressure button (air pressure start/stop)

- After turning on the machine, press the air pressure switch button to turn on the air pressure function. The corresponding air pressure indicator light is bright green, the default air pressure is low, and the corresponding air pressure intensity light is green.
- Press the button again to turn off the air pressure function. The corresponding air pressure LED light goes out and the air pressure intensity light goes out.
- Air pressure intensity button (low/medium/high)After turning on the air pressure function, press this button to turn on the low air pressure, and the corresponding LED light will turn green.
- Press the button again to turn on the medium air pressure, and the corresponding LED light will turn yellow.
- Press the button again to turn on the high air pressure, and the corresponding LED light will turn red.
- This cycle repeats.
- The air pressure intensity button can only be enabled when the air pressure button is turned on.

Heating button (low/high)

- After turning on the machine, press this button to perform low-speed heating, and the corresponding green LED light will light up. Press this button again to perform high-speed heating, and the corresponding red LED light will light up.
- Press this button again to turn off heating, and the corresponding red and green lights will go out.

Mode button (A/B/C mode)

- After turning on the power, press this button to execute automatic mode A, and press the button again to execute automatic mode B. Press the button again to execute automatic mode C. Press the button again to execute automatic mode A, and so on.

During automatic mode, you can operate the heating and gas pressure functions through the corresponding buttons. If the heating is not turned on before the mode is turned on, pressing this button will turn on the low heating by default (green light).

- If the gas pressure is not turned on before the mode is turned on, pressing this button will turn on the medium gas pressure by default (yellow light).

ATTENTION

Please use a low-speed massage function when you start using it. If you feel any abnormality during use, please stop using it. It is recommended that you do not use it for more than 30 minutes each time.

Product Specifications

Model	EMK-795A
Name	Foot Massager
Power input	DC 24.0V \approx 2A
Power consumption	48W
Default Time	15 mins

This product is intended for house-hold use only.

ATTENTION

Note: Certified manufacturer's adapter must be equipped in the use of this products.

How to clean

Storage method

- To keep the product clean, make sure it is away from highly humid places. You should store the device in its original packaging and keep it in a clean, dry place.
- Do not place it in direct sunlight or high temperature.
- If it is not used for a long time. Please cover it with a dust cover and unplug the power plug from the power socket.

1. Cleaning

Before you clean, please make sure the power is switched off, and the power plug has been pulled out from the socket.

Please use dry cloth to clean the product.

2. Maintenance


Keep the product clean and far from the wet environment.

Clean when you feel cleaning is necessary.

Cleaning the main unit

Wipe away dirt with a cloth soaked in soapy water, wring out the cloth before wiping.

Do not use alcohol, paint thinner, benzene, or nail polish remover. Using these may damage the appliance or cause cracks or discoloration in the components.

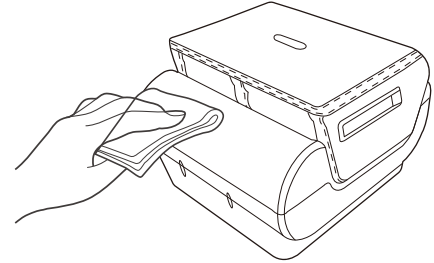
Q&A

Problem	Cause	Remedy
The sound of the product working normally or the sound of "whooshing" or "tapping" occurs during the operation.	The sound made when the product is working normally or the sound of friction between the massage ball and the cloth cover.	This does not indicate a problem.
Unable to turn on the device.	Set time has expired.	Use the device after 30 minutes of shutting down.

If you still have any problem after using the methods above, please contact our customer support team.

Correct Disposal of this product

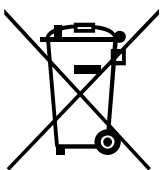This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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