IN240701680V01

712-067V00

EN_ IMPORTANT, RETAIN FOR FUTURE REFERENCE: READ CAREFULLY.

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- **INTRODUCTION**

The rollators are manufactured from high-quality materials and include a seat, bag, and locking brakes. The hand-bar are adjustable in height and the rollator folds easily for transport or storage.

- **INTENDED USE**

Designed to assist the user to walk, rest and carry small items of shopping for both indoor and outdoor use on flat surfaces by a single individual with limited mobility with a maximum weight as shown. The intended user should possess the necessary physical, visual, and cognitive skills to operate the rollator safely. If you have any doubt, please contact qualified technician. It is important that the user reads and understands the safety advice given within this manual before attempting to operate the product to avoid possible injury to the user or bystanders and avoid damage to the rollator or surrounding property.

- **SAFETY ADVICE**

This rollator is designed to be used as a mobility aid, creating a natural gate as the user walks securely with their primary weight between the handle bars with the unit in front of the user. Read this manual carefully before using your rollator.

Your rollator is equipped with various components, you should know these components and parts before proceeding to read this user manual. Design and specifications could be changed without further notices.

• SPECIFICATION

Overall width 65.5 cm

Overall depth 75cm Seat area 46*25 cm

Adjustable handle height 80-92.5cm

Max. user weight 136 kg

Net weight 8.3 kg

Assembly instructions

1. Remove the rollator from the packing and remove the protective foils. Spread out the parts for assembly as shown.

2. Take out the main frame and press the both sides of the seat, then the seat was secured at a click.

3. Remove the rear wheel with marking "RB", and insert the rear wheel with marking "RB" into the frame with marking "RB"

4. Remove the rear wheel with marking "LB", and insert the rear wheel with marking "LB" into the frame with marking "LB"

5. Remove the wing screw on the front wheel

6. Insert the front wheel into the frame

7. Fasten the front wheel with the previously loosened wing screw. The Black wing screw pointing outwards. When turning, the screw head will tighten in the frame. And then, Repeat steps 5.- 7. for these cond front wheel on the other side of the frame.

8. Our recommendation: Fix the brake cables on both sides of the frame using. The locking buckle as shown.

9. Attach the bag. Hook the tabs of the bag into the Hooks provided.

Adjusting the height

Press the button on the height adjustment on the frame of the rollator and adjust the handle to the appropriate height. Repeat for the other side. Make sure the buttons on both handles are locked. Make sure that both handles are set to the same height.

Handle height: 80-92.5 cm

Your ideal handle height

Adjust the handles of the rollator so that they are at the height of your wrist when your arms are hanging. When you place your hand on the handle, your elbow should be slightly bent. If you find that you need to bend over or keep your arms bent at a right angle, adjust the height of the handles.

Brake

1. Pull the brake handles up to brake the rollator.

2. Push the brake handles down to lock the rear wheels. In this position, the parking brake of the rollator is active. To release the parking brake again, pull the brake handles upward.

3. Brakes must be applied in the lock position when the seat is going to be used. Precheck the lock by pushing forward on the rollator.

Brake adjustment is carried out as follows:

Turn the adjusting screw anti-clockwise to increase the tension, then lock the adjustment locking nut clockwise whilst holding the adjuster screw stationary. Brake adjustment reduced to allow light pressure may be ineffective on the brake position.

Removable cane holder

1. Attach the cane holder to the side of the frame in the provided position.

2. Attach the cane holder to the side of the frame in the provided position. To do this, slide the cane holder onto the holder from above.


Folding and storing the rollator

1. Pull the handle in the middle of the seat firmly upwards to fold the rollator.

2. Continue to push the rollator on both sides of the frame.

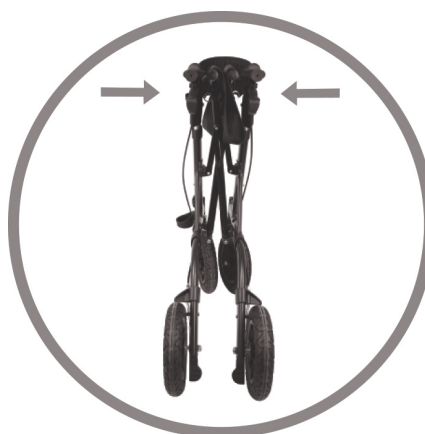MAINTENANCE

Over a period of time, the brake cables may stretch, and if no further adjustments can be made, the cables should be replaced. Seek the help of the manufacture when replacing brake cables. The castors should be checked that they are secure each time before using the rollator. Also, check the brakes operate correctly and the handles and hand knobs are secure. Should any part of the rollator appear loose or insecure, discontinue use until it is inspected by a qualified technician. The following shows the maintenance schedule:

Weekly

- Check brake for malfunction
- Check all wheels and castors for loose bolts and nuts and knobs

Periodically

- Check hand grips for wear
- Check frame bolts are secure
- Check front and rear axle bolts are secure
- Check tires for wear or cracks

Yearly

Once a year have your rollator fully inspected, cleaned, maintained and lubricated with all moving parts

Serviceable Components

Wheels

- Hand grips
- Bag
- Knobs
- Cane holder
- Standing armrest sleeve

CLEANING

Your rollator can be cleaned with a damp cloth or sponge and warm soapy water. Domestic anti-bacterial surface cleaners may be used if required. Care should be taken not to leave any cleaning residues on the hand grips or seat(wherfitted)which may leave them slippery when used. Do not use as couring pad for scrapers this may damage surface coatings.

WARNINGS!

- 1.Maximum working limit 136kg.
- 2.Do NOT stand on this product.
- 3.Ensure both push handles are adjusted to the same height.
- 4.Seek advice from a trained healthcare professional for height adjustment of this product.
- 5.Always maintain a safe distance between your rollator and other pedestrians or objects.
- 6.Some components may cause harm or damage when they come into contact with other objects.
- 7.Do not attempt to replace the wheels or castors of your rollator. If they are faulty or require a replacement, please contact qualified technician.
- 8.Ensure that your rollat or is opened fully(folding scissor armis straight)before reusing it.
- 9.Ensure that your fingers are clear from the scissor arm when folding.

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If you have any questions, please contact our customer care center, Our contact details are below:

US

 001-877-644-9366

 customerservice@aosom.com

Imported by Aosom LLC
7777 Center Ave. Suite 430
Huntington Beach, CA 92647-9998 USA
MADE IN CHINA

If you have any questions, please contact our customer care center, Our contact details are below:

CA

 416-792-6088

 customerservice@aosom.ca

Imported by Aosom Canada Inc.
Unit 33-1 Sparks Avenue, North York, Ontario,
Canada
M2H 2W1
MADE IN CHINA

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If you have any questions, please contact our customer care center, Our contact details are below:

UK

 0044-800-240-4004

 enquiries@mhstar.co.uk

IMPORTER ADDRESS:
MH STAR UK LTD
1 Northampton Cross Logistics Park
Northampton Cross Road
Northampton
Nn4 9FH
Tel: 0800 240 4050
enquiries@mhstar.co.uk
MADE IN CHINA

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