**UK
CA** **CE**

9
STEPS

0.5 HR

EN_ IMPORTANT, RETAIN FOR FUTURE REFERENCE: READ CAREFULLY.

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
1. Components of the rollator

In order to read and understand these instructions, it is important that you familiarize yourself with the names of the various components of the rollator. The figure below shows you an overview of your rollator and its components.

- 1. Backrest
- 2. Brake line
- 3. Push handle
- 4. Adjustable parking/running brake
- 5. Nylon buckle strap for attaching the walking stick
- 6. Red button of handle height adjustment
- 7. Screw for adjusting the brake
- 8. Rear wheel
- 9. Front wheel
- 10. Bag
- 11. Nylon buckle strap for securing the brake line
- 12. Cane holder
- 13. Foldable seat

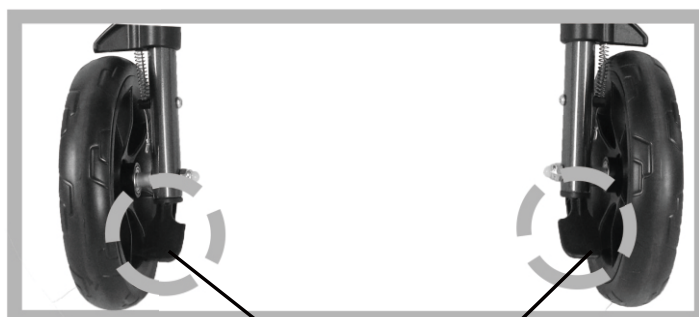Step aid

2. Scope of Delivery

1. Left and right rear wheels with left and right handle with backrest
2. Two front wheels
3. Main frame with seat and Two Nylon buckle straps for securing the brake lines
4. Bag
5. Cane holder


3. Assembly of the rollator

- 3.1. Remove the rollator from the packing and remove the protective foils. Spread out the parts for assembly as shown.

- 3.2. Take out the main frame and press the both sides of the seat, then the seat was secured at a click.

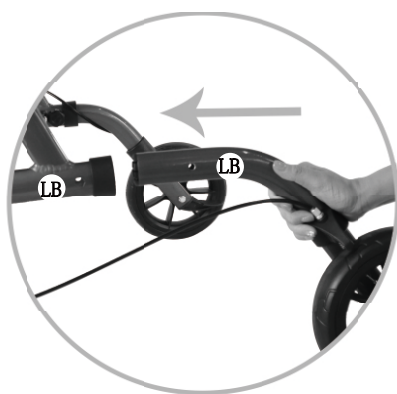- 3.3. Remove the Plum-shaped nuts and bolts on the rear wheel. Take the left rear wheel with marking "LB" insert into the frame with marking "LB" and take the right rear wheel with marking "RB" insert into the frame with marking "RB". Fasten the rear wheels with the previously loosened Plum-shaped nuts and bolts.

3.4 Remove the Plum-shaped knob on the front wheel.
Insert the front wheel into the rollator frame from the front.

3.5. Fasten the front wheel with the previously loosened Plum-shaped knob.
The Plum-shaped knob is pointing inward. When turning, the knob will tighten in the frame. And then, Repeat steps 3.4. - 3.5. for the second front wheel on the other side of the frame.

3.6. Press the red button and insert the handle into the rollator frame.
Make sure that the two handles are at the same desired height.

3.7. Our recommendation: Fix the Nylon buckle on both sides of the frame using. The Nylon buckle as shown.

3.8. Attach the bag. Hook the bag into the hooks provided.

3.9. Fastening the cane holder.
Slide the cane holder from above onto the holder on the left of the rollator frame.

4. Handling/Operating

In order to use the rollator safely, you need to make some adjustments and check them regularly.

4.1. Adjusting the height

Press the red button of the height adjustment on the frame of the rollator and adjust the handle to the appropriate height.

Repeat for the other side. Make sure the red buttons on both handles are locked. Make sure that both handles are set to the same height.

To determine the correct height of the handles, stand upright and let your arms hang loosely. In this position you must be able to loosely grasp the handles.

Only use the rollator with the handle height correctly adjusted for your height.

Only with a correctly adjusted handle height you can avoid incorrect loads and injuries.

4.2. Operating the brakes

Always operate both handbrake levers at the same time, as each brake only acts on one wheel. Be careful not to crush your fingers when applying the brakes.

Before each use, check that the brakes are working properly.

To apply the parking brake, press down the two brake handles and engage them. Check that the rollator is stationary before sitting down.

To release the brakes again, guide the two brake handles upwards.


Important!

If you want to brake while walking, pull the brake lever up slightly to reduce your speed.

Do not push the rollator when the brakes are firmly applied or locked.

This can cause damage to the brakes or the wheels.

Failure of the brake system may result. Do not use the rollator if the brake system is defective.

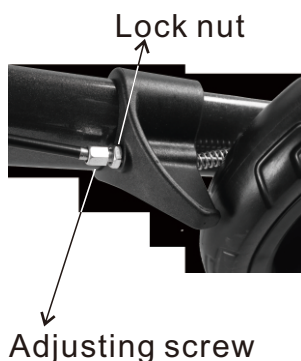4.3. Brake adjustment is carried out as follows:


To adjust the rear brake, adjust the adjusting screw and locking nut. These are located on the rear wheels of the rollator. Turn the adjusting screw counterclockwise to increase the tension. Then secure the adjustment by turning the locking nut clockwise. Hold the adjusting screw in place while doing so.

If the rear brakes are set too lightly or too tightly, the braking effect may be impaired or the wheels or brakes may be damaged.

4.4. Folding and storing the rollator

4.4.1. Pull the strap in the middle of the seat firmly upwards to fold the rollator.

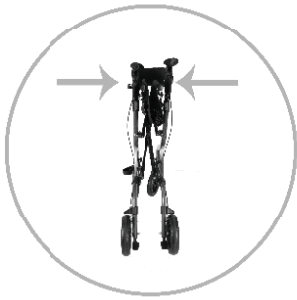4.4.2. Continue to push the rollator on both sides of the frame.

4.5. Overcome obstacles

- Hold the rollator firmly when negotiating obstacles.
- Take special care when walking down kerbs.
- Do not use the rollator on stairs.
- Do not tilt the rollator more than 5 degrees to one side.

4.6. Step aid

If you have to overcome kerbs or similar obstacles, step on the step aid and pull the rollator towards you by both handles to lift the front wheels. Approach obstacles at an angle and lift the front wheels over the obstacle one after the other.

4.7. Using the seat correctly

- Only sit on the rollator when both parking brakes are applied.
- Sit in the middle of the seat.
- Make sure that both feet are touching the ground.
- Do not sit on the rollator if the ground is very slippery. Even though the parking brakes are applied, the rollator may slip away.
- The seat may only be used by one person at a time.

4.8. Walking correctly with the rollator

- The rollator must always be held close to the body when walking so that it can support you well and enables a safe gait.
- If the rollator is too far away from the body when walking or if you walk too bent over, the rollator can easily roll away to the front.

Danger of falling!

- Always walk upright and keep your gaze forward.
- Only use the rollator if the handle height is correctly adjusted for your height (see chapter "Adjusting the handle height").

5. MAINTENANCE

Check the functionality of the rollator before each use.

To do this, check the following points:

- The brakes are in perfect working order.
- The brake cables are free of deformations and cracks.
- The wheels are in perfect working order. All wheels can be turned properly.
- The wheels do not show any heavy signs of wear.
- All screw connections are tight.
- The entire rollator is free of damage, cracks and coarse dirt.
- The folding mechanism of the rollator under the seat is free of damage and dirt.
- The folding mechanism works perfectly.

6. Safety notice

- Before each use of the rollator, check that the brakes are working properly.
- Make sure that the rollator is properly engaged when you unfold it.
- Before each use of the rollator, check that the push handles are properly locked.
- Always apply the parking brake when you take a break or park the rollator.
- Make sure that the correct push handle height has been set to avoid incorrect loads and falls. Do not push the rollator too far in front of you.
- Watch out for crushing and shearing points between the moving parts of the rollator. Be especially careful when folding or unfolding the rollator and when adjusting the height of the push handles.
- Walk with the rollator on the pavement and not on the roadway. Make sure that you are seen. Wear light-coloured clothing and use reflectors when using the rollator outdoors.
- Only use the rollator for the purpose for which it was designed. Do not use the rollator to transport loads.
- Pay particular attention to children in the vicinity of the rollator. Do not allow children to play with the rollator.
- The rollator may only be used by persons with a maximum weight of 136kg.
- The bag may be loaded with a maximum of 5kg.
- The rollator must be checked regularly. All screw connections and joints must be tight. Components must not be damaged or loose.
- In case of damage or malfunction, further use of the rollator is prohibited.
- Do not make any changes or modifications to the rollator. Changes or modifications can negatively affect the safety of the rollator.
- This product should avoid prolonged exposure to direct sunlight or high-temperature environments to prevent aging, fading, performance degradation, or component damage.
- For storage, please keep it in a cool and dry place.

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If you have any questions, please contact our customer care center, Our contact details are below:

US

 001-877-644-9366

 customerservice@aosom.com

Imported by Aosom LLC
7777 Center Ave. Suite 430
Huntington Beach, CA 92647-9998 USA
MADE IN CHINA

If you have any questions, please contact our customer care center, Our contact details are below:

UK

 0044-800-240-4004

 enquiries@mhstar.co.uk

IMPORTER ADDRESS:
MH STAR UK LTD
1 Northampton Cross Logistics Park
Northampton Cross Road
Northampton
Nn4 9FH
Tel: 0800 240 4050
enquiries@mhstar.co.uk
MADE IN CHINA

If you have any questions, please contact our customer care center, Our contact details are below:

CA

 416-792-6088

 customerservice@aosom.ca

Imported by Aosom Canada Inc.
Unit 33-1 Sparks Avenue, North York, Ontario,
Canada
M2H 2W1
MADE IN CHINA

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