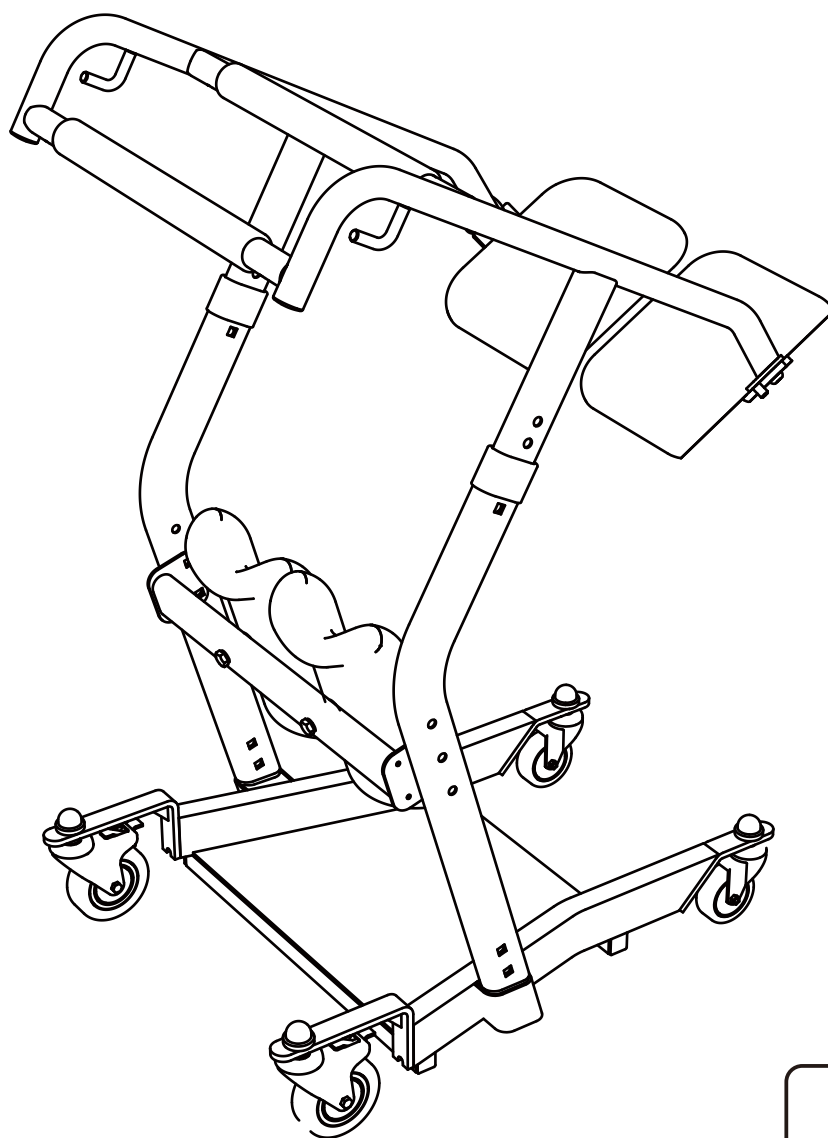
5
STEPS

0.5 HR

EN_ IMPORTANT, RETAIN FOR FUTURE REFERENCE: READ CAREFULLY.

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






Thank you for choosing our Sit to Stand Patient Lift !



Product Structure Diagram

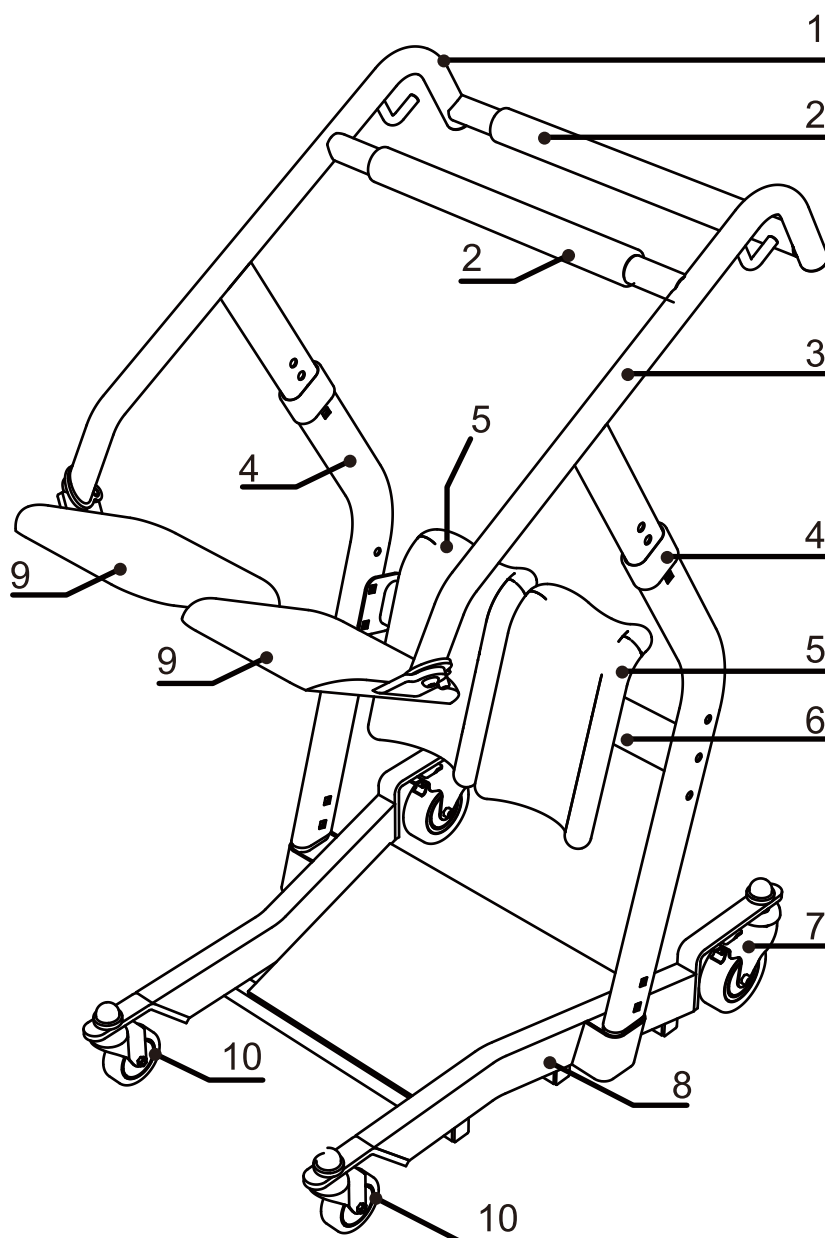- 1-Left armrest assembly
- 2-Handrail tube
- 3-Right armrest assembly
- 4-Oval hollow sleeve
- 5-PU leg pads
- 6-Leg pad fixing assembly
- 7-Front wheel assembly
- 8-Base frame assembly
- 9-PU seat cushion
- 10-Rear wheel assembly

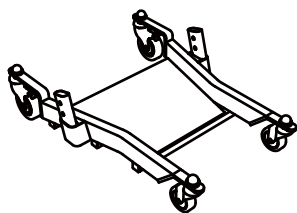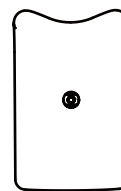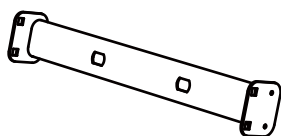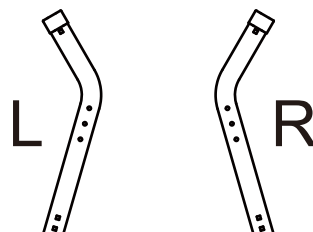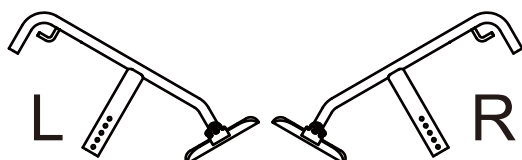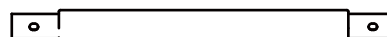A x 1**B** x 2**C** x 1**D** x 1**E** x 1**F** x 2

HARDWARE

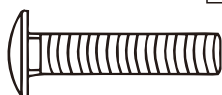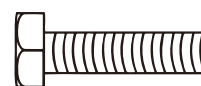1 x10
M10*45**2** x12**3** x10**4** x10**5** x 2**6** x 2
M10*25**7** x 2**8** x 4
M6*15**9** x 1**10** x 1

Preparations before assembly

After opening the packaging box, first use scissors to remove the plastic ties and bubble wrap. Please take care to protect all parts and components inside the box. If assembling on the ground, first lay down cardboard or carpet to prevent scratches. Be sure to distinguish between parts and components, then remove the inner cardboard box and arrange the parts neatly. The packaging box contains loose parts as shown in the Figure above.

Installation Instructions

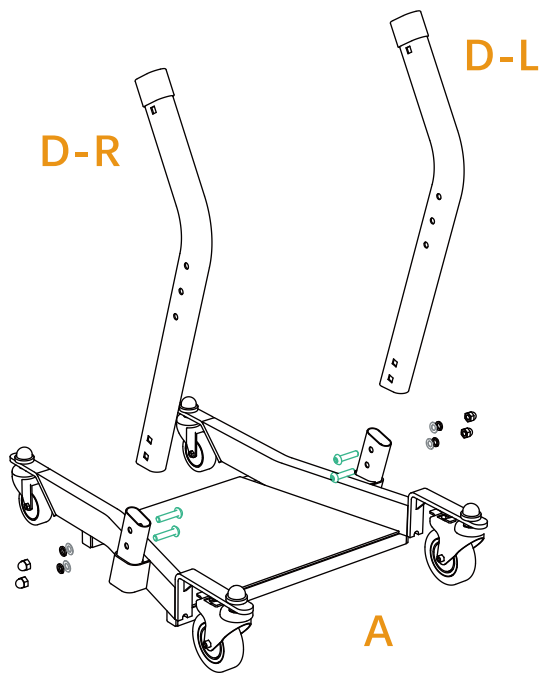
Figure 1

STEP 1

As shown in Figure 1, align the D-L left support elbow and D-R right support elbow with the hole in the A. base frame assembly, insert the pipe, and secure it with M10*45 carriage bolts and 10.2 washers. Then use an open-end wrench to tighten the M10 cap nut.

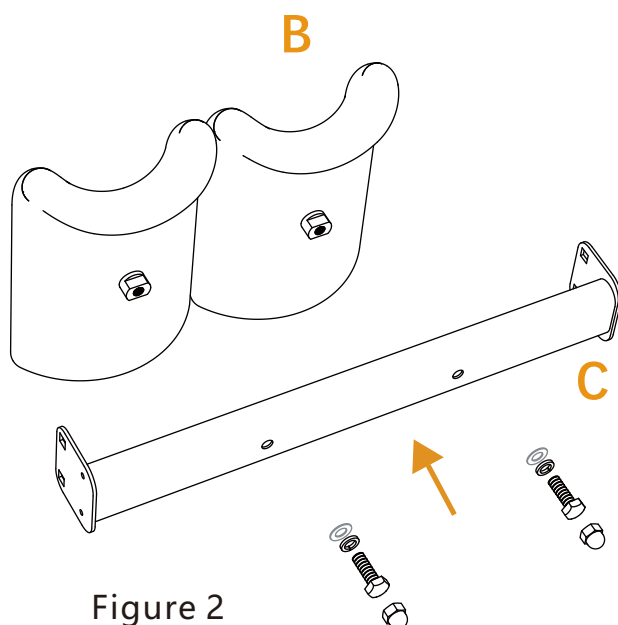
Figure 2

STEP 2

As shown in Figure 2, use accessory B leg pads and C leg pad fixing components to insert M10 carriage bolts into the holes, then place 10.2 spring washers and 10.2 flat washers, and tighten with an open-end wrench and M10 cap nuts.

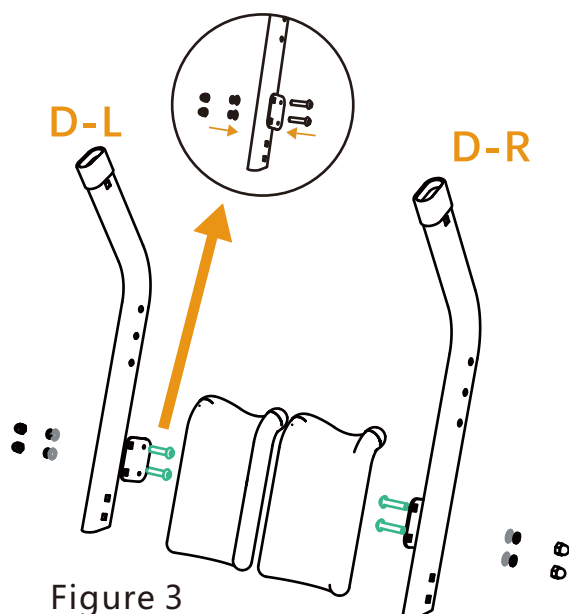
Figure 3

STEP 3

As shown in Figure 3, align the assembled C. leg pad fixing assembly with the holes on the D-L left support elbow and D-R right support elbow pipes, insert the M10 carriage bolts into the holes, and then place the 10.2 spring washers and 10.2 flat washers. Use an open-end wrench to tighten the M10 cap nuts.

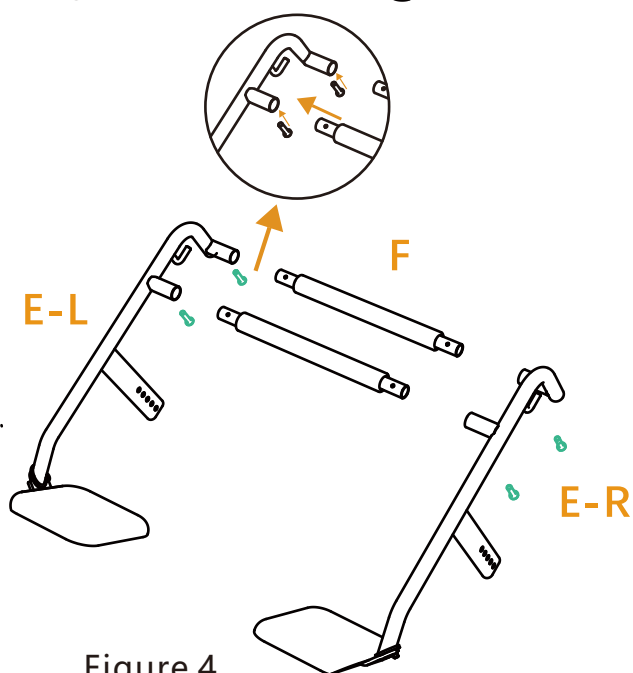
Figure 4

STEP 4

As shown in Figure 4, insert the F handrail tube into the E-L left handrail frame and E-R right hand-rail frame tubes, then align the holes in the tubes and use a hex wrench to tighten the M6*15.

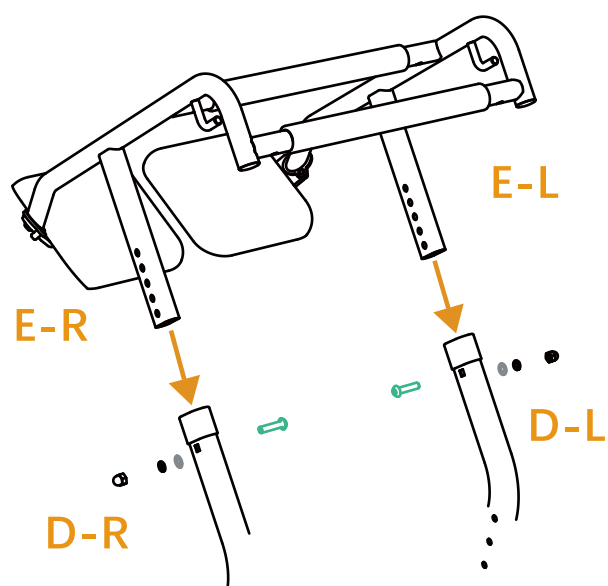
Figure 5


STEP 5

As shown in Figure 5, insert the E-Left armrest frame and E-R right armrest frame into the aligned holes inside the D-L left support elbow pipe and D-R right support elbow pipe. Insert the M10 carriage bolts into the holes, followed by the 10.2 spring washers and 10.2 flat washers. Use an open-end wrench to tighten the M10 cap nuts.

Post-Installation Inspection

1. Refer to the instructions in the manual to check whether the standing frame is installed correctly.
2. Verify that all accessories for the standing frame are assembled completely, and check for any extra accessories.
3. Check whether the standing frame is loose. If so, tighten the screws and nuts.
4. Inspect all components of the standing frame for any misalignment and make appropriate adjustments.
5. Check whether the standing frame brakes are working properly, whether they fold and unfold smoothly, and whether they turn smoothly.

Corresponding Height

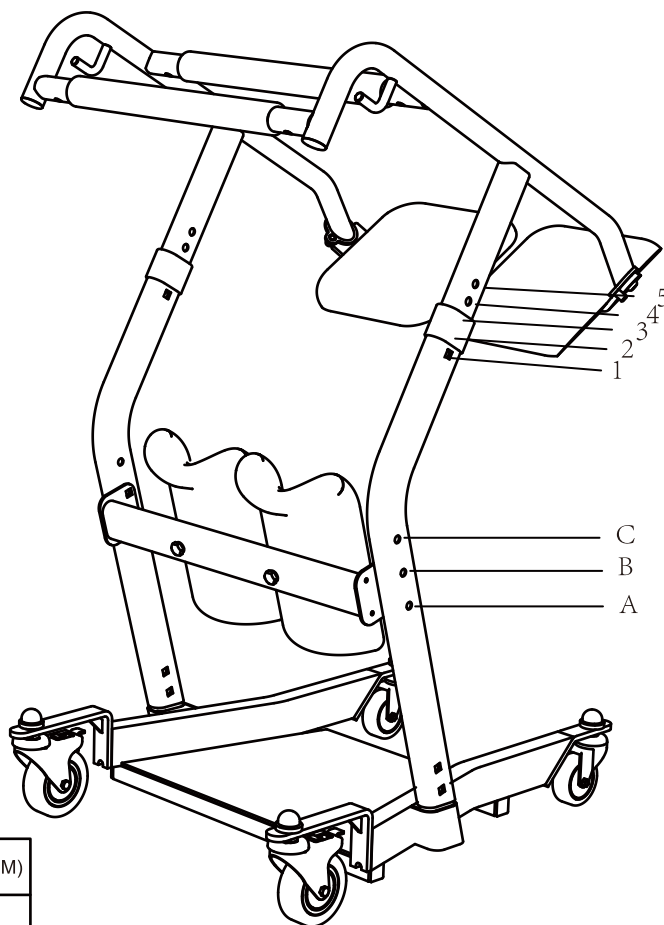
Leg pad adjustment	Seat cushion adjustment	Corresponding Height(CM)
AB	1	150-160
	2	155-165
	3	160-170
	4	165-175
	5	170-180
BC	1	155-165
	2	165-175
	3	170-180
	4	175-185
	5	180-190


Please note that this height correspondence is not absolute.
For specific adjustments, please refer to the actual situation.

How to Use


1. When using the standing frame, the operator first pushes the standing frame to the front of the user's bed, steps on the brakes on both sides of the front wheels, and the user stands on the non-slip foot pads on the base frame of the standing frame, holding the handrails with both hands, as shown in the figure below:

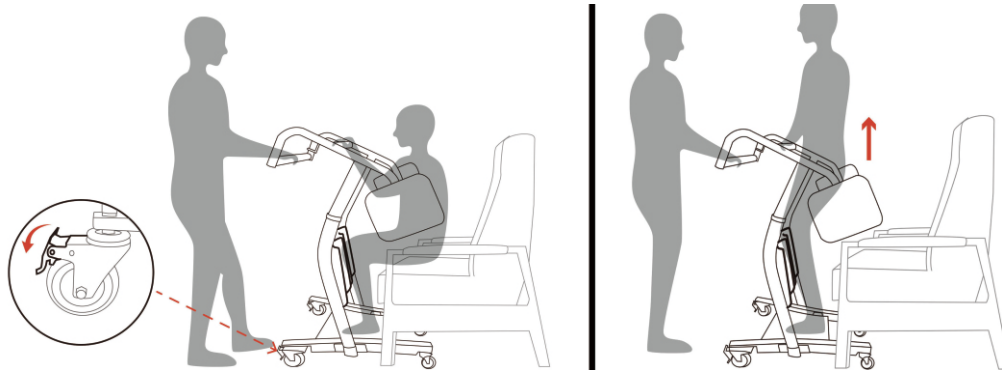2. Turn the upper PU seat cushions of the two armrests inward toward the user's back. The user sits down, holds onto the front armrests, and then releases the two front wheel brakes. The operator can then pull the user to move the frame, as shown in the figure below:

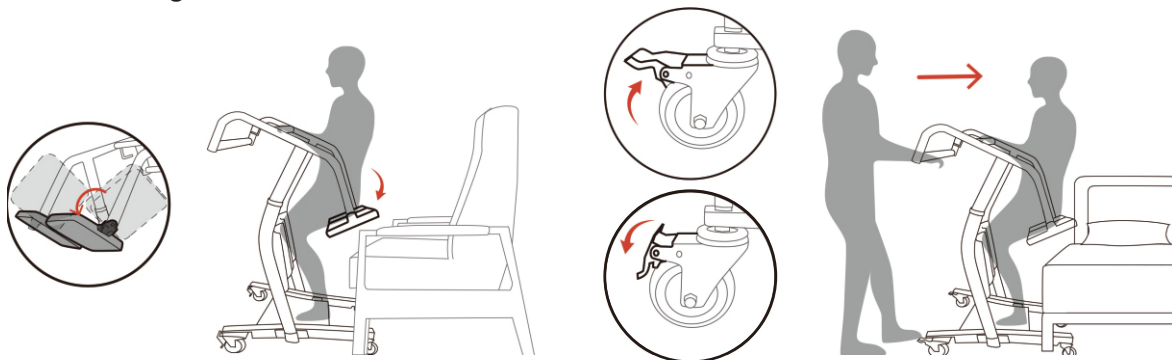3. When the operator needs to transfer the user to the bed, push the user to the front of the bed, apply the brake, have the user stand up, flip the PU seat cushion of the handrail frame outward on both sides, and slowly remove the standing frame after the user is seated securely on the bed, as shown in the figure below.

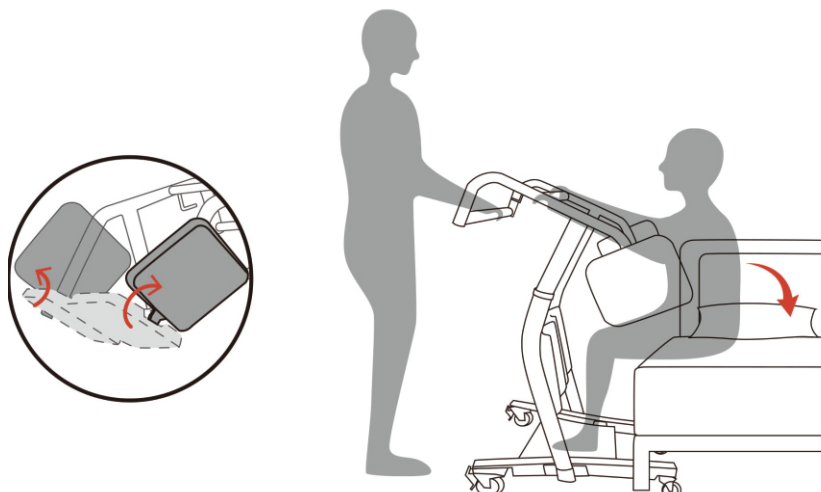Product Maintenance and Care Methods

1. Please do not scratch or rub the product surface with hard objects to prevent damage to the material.
2. For daily cleaning and maintenance, use a neutral synthetic detergent, wipe with a damp cloth, and then wipe dry with a dry cloth.
3. Avoid contact with strong acids and alkalis to prevent corrosion.

Storage Conditions and Methods

After packaging, the standing rack should be stored at a temperature not exceeding - 10°C to 40°C, with relative humidity not exceeding 80%. It should be stored in a clean indoor environment with good ventilation and no corrosive gases. During transportation, it should be protected from moisture, direct

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If you have any questions, please contact our customer care center, Our contact details are below:

US

 001-877-644-9366

 customerservice@aosom.com

Imported by Aosom LLC
7777 Center Ave. Suite 430
Huntington Beach, CA 92647-9998 USA
MADE IN CHINA

If you have any questions, please contact our customer care center, Our contact details are below:

UK

 0044-800-240-4004

 enquiries@mhstar.co.uk

IMPORTER ADDRESS:
MH STAR UK LTD
1 Northampton Cross Logistics Park
Northampton Cross Road
Northampton
Nn4 9FH
Tel: 0800 240 4050
enquiries@mhstar.co.uk
MADE IN CHINA

If you have any questions, please contact our customer care center, Our contact details are below:

CA

 416-792-6088

 customerservice@aosom.ca

Imported by Aosom Canada Inc.
Unit 33-1 Sparks Avenue, North York, Ontario,
Canada
M2H 2W1
MADE IN CHINA

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